

THROMBO inCode


More answers for managing
genetic thrombophilia

GEN inCode

HEREDITARY THROMBOPHILIA

Hereditary thrombophilia is a genetic predisposition to venous thromboembolism (VTE).

It is due to several genetic¹ mutations that alter the clotting cascade, increasing individual risk² of VTE.

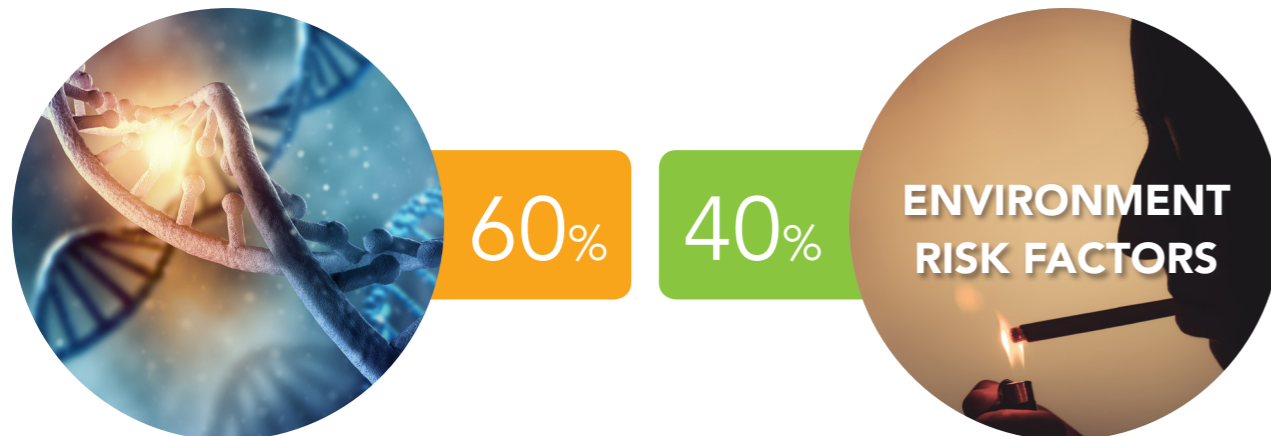

According to the GAIT study, genetics accounts for 60% in the occurrence of thromboembolic events.³

Genetics, along with plasma tests, is key to diagnosing thrombophilia.

Unlike plasma tests, the results of genetic tests are not affected by clinical or pharmaceutical conditions and they don't change over the course of a patient's life.

CURRENT DIAGNOSIS OF HEREDITARY THROMBOPHILIA⁴

PLASMA TESTS	+	GENETIC TESTING
<ul style="list-style-type: none"> • Prothrombin time • Activated partial thromboplastin time • Antithrombin activity • Protein C activity (chromogenic) • Free protein S antigen 		<ul style="list-style-type: none"> • Factor V Leiden • 20210 AG PT gene mutation (Prothrombin)

The genetic markers analysed today, FV Leiden and Prothrombin, can only diagnose 20% of patients with hereditary thrombophilia.⁵

MORE COMPLETE RESPONSE TO MANAGING THROMBOPHILIA

THROMBO inCode combines:


- 1 A validated panel⁵ of 12 genetic mutations, all of which are causal and functional with a pro-thrombogenic effect on the clotting cascade.

Screening for thrombophilia genetic profile.

- 2 An algorithm using the patient's genetic profile and medical history/risk factors: gender, age, BMI, diabetes, smoking habits, family history of VTE, hormone-based contraceptives, pregnancy. The different risk levels are calculated based on various combinations of potentially modifiable risk factors.

Calculation of the patient's risk of having a venous thromboembolic event.

THE TOOL FOR GENETIC DIAGNOSIS OF HEREDITARY THROMBOPHILIA AND PREDICTION OF VTE RISK

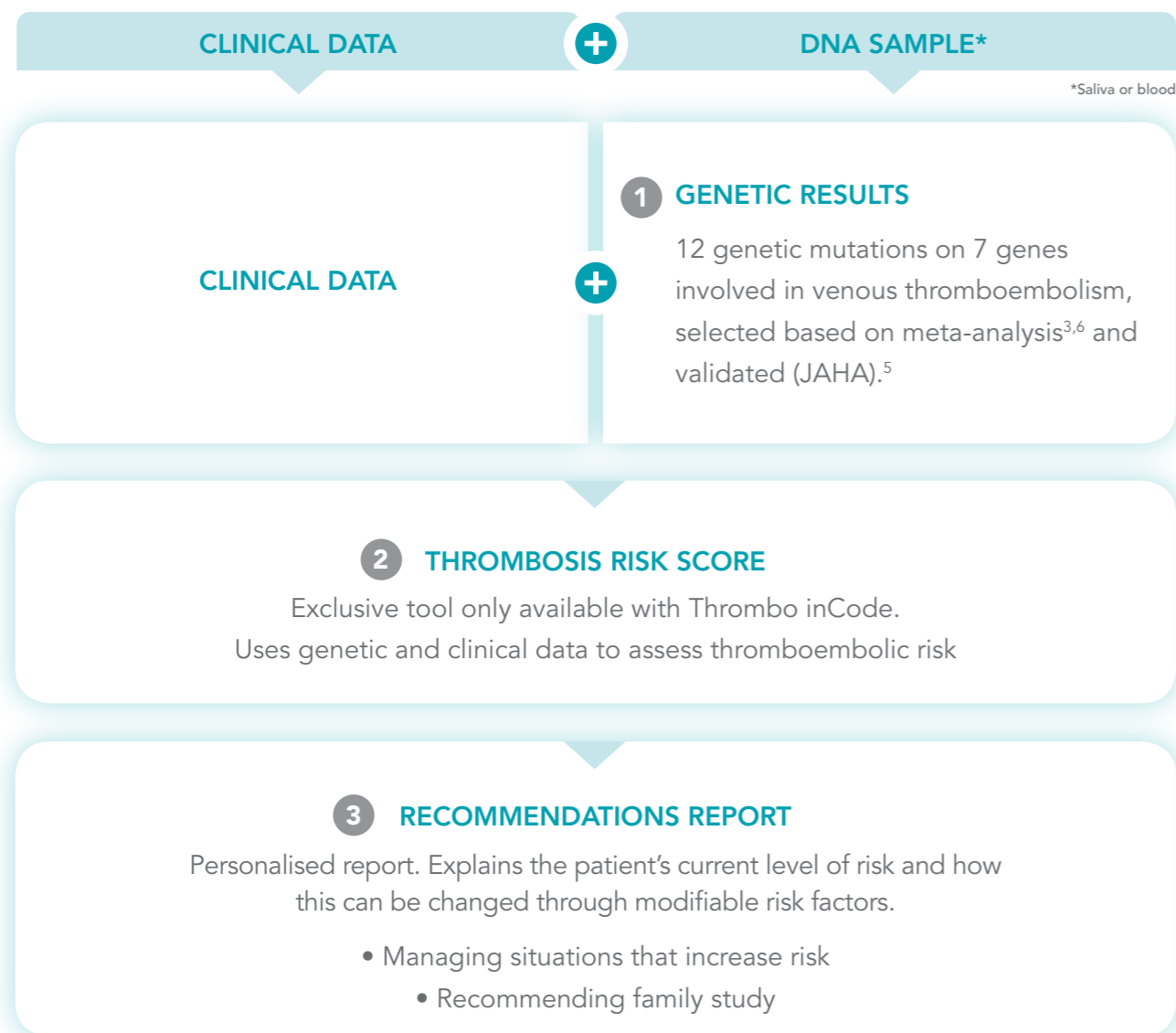


THROMBO inCode
IS AVAILABLE AS A DIAGNOSTIC SERVICE OR LABORATORY KIT

THE PANEL OF 12 GENETIC MUTATIONS

VARIANTS	GENE/PROTEIN MUTATIONS	PREVALENCE	RELATIVE VTE RISK ^d
rs6025 (FV Leiden)	Factor V ⁷⁻¹³	15-25%	5 (h) / 34 (H)
rs1799963 (Prothrombin)	Factor II ^{7-9, 14,15}	6-16%	4,20 (h) / 26,4 (H)
rs1801020	Factor XII ¹⁶⁻¹⁸	6% ^c	3,22 (h)
rs8176719	ABO Group ¹⁹⁻²¹ (A1 carriers)	nd	2,6
rs7853989			
rs8176743			
rs8176750			
rs2232698	Serpin A10 ²²	4,4%	3,89
rs121909548	Serpin C1 ²³	1,7%	9,75
rs118203906 (FV Cambridge)	Factor V ⁷⁻¹³	nd	5 (h) / 34 (H)
rs118203905 (FV Hong Kong)			
rs5985	Factor XIII ^{8,24,25}	2% ^c	0,8

(h): Heterozygote / (H): Homozygote / c: Zygosity / d: Data obtained from bibliography and meta-analyses
VTE: Venous thromboembolism

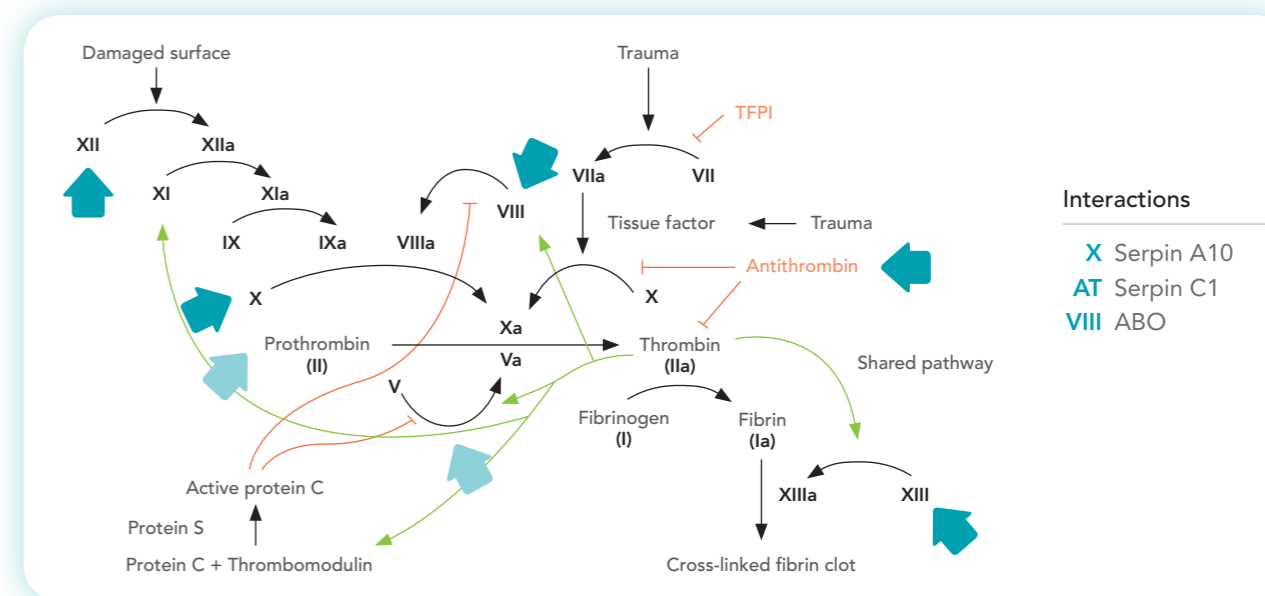



VALIDATION STUDY⁵

Published in the JAHA (Journal of the American Heart Association) • October 2014

AIM To analyse the predictive capacity of ThromboinCode to assess a subject's risk of having a thromboembolic event.

STUDY POPULATION

SANT PAU – Spain

- 248 cases – 249 controls
- Age: 49.0/47.1
- Males: 44.0/44.6%
- Representative sample of the general population (Caucasian)

MARTHA – France

- 477 cases – 477 controls
- Age: 44.2/43.9
- Males: 41.5/41.5%
- Sample population (Caucasian) enriched with PT and FVL

Comparison of discriminative and predictive capacity of the 12 mutations in Thrombo inCode (GRS2) with:

- GRS1** V Leiden factor and 20210 AG Prothrombin gene mutation (current gold standard for diagnosing hereditary thrombophilia)
- GRS3** Risk score developed by Dr Rosendaal's thrombophilia research group (Leiden University)
- GRS4** Combination of GRS2 and GRS3, taking into account genetic mutations associated with the development of venous thromboembolism, according to the literature.


*GRS: Genetic Risk Score

VALIDATION STUDY⁵ RESULTS

THROMBO inCode

- ▶ **SIGNIFICANTLY IMPROVES** the ability to predict a VTE event based on family history (AUC 0.701 vs 0.589).
- ▶ **SURPASSES** FVL + PT combination in ability to predict a VTE event.(AUC 0.68 vs 0.57) and increases sensitivity from 20% to 85%.
- ▶ Adding **more variables does not increase** the predictive capacity of Thrombo inCode.

Clinical utility (measured in terms of sensitivity) of TiC* compared to FVL + PT in SANT PAU population

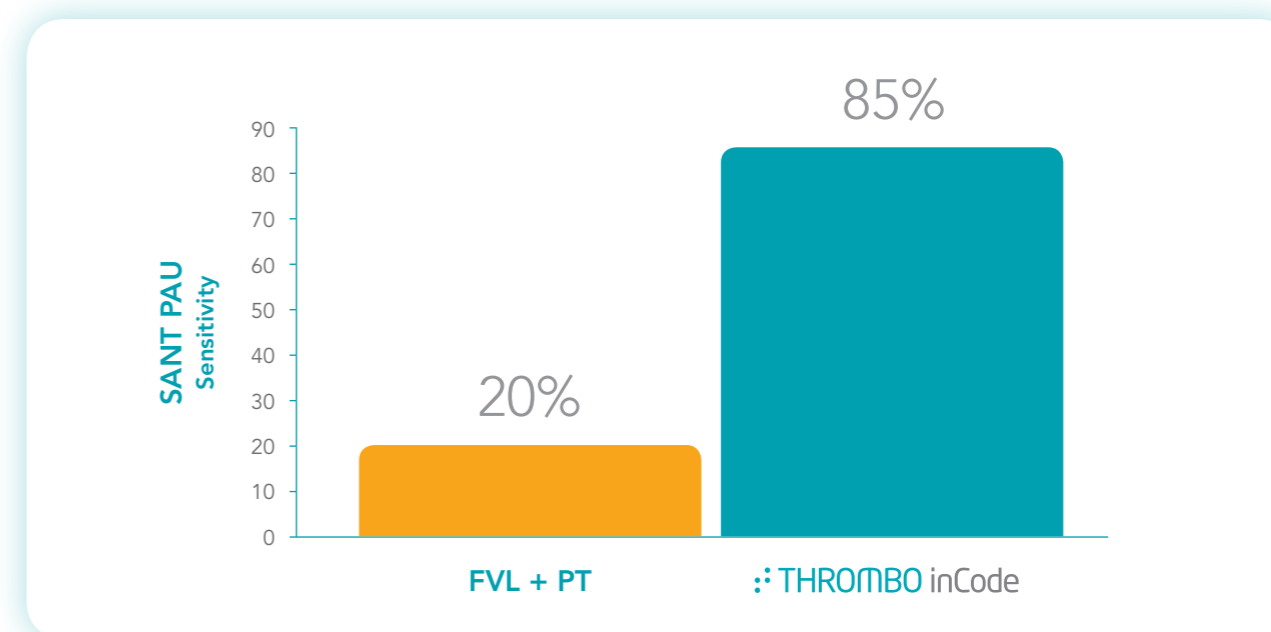

*TiC: Thrombo inCode / Graph adapted from reference text 10

ALGORITHM TO CALCULATE VENOUS THROMBOEMBOLISM RISK SCORE

This algorithm allows us to calculate a person's probability of having VTE depending on their genetic, sociodemographic and clinical characteristics.



- It weighs each value taking into account its associative power (OR*). Each one has its own specific weight.
- Most algorithms used today only take into account clinical data, while **Thrombo inCode uses both genetic and clinical data.**
- The risk assessment also takes into account specific conditions that may change over time, allowing it to be recalculated whenever necessary.

*OR: Odds ratio

$$X_1 \beta_1 + X_2 \beta_2 + X_3 \beta_3 \dots + X_n \beta_n$$


VTE RISK SCORE

Each variable (x1, x2, x3, etc.) in the algorithm corresponds to clinical or genetic conditions assessed in our algorithm:

- Genetic mutations (PT, FVL, FXII, FXIII, etc.)
- Age, gender, smoking habits, diabetes, etc.

Each variable is given its own specific weight ($\beta_1, \beta_2, \beta_3$, etc.). This ensures it not only takes into account the presence of the risk factor but also quantifies its degree of risk.

ALGORITHM TO CALCULATE VENOUS THROMBOEMBOLISM RISK SCORE

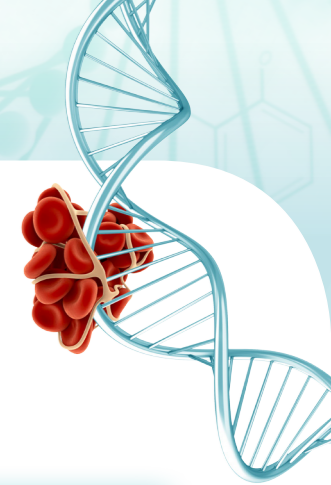

FICTITIOUS PATIENT					
RISK SCORE BASED ONLY ON CLINICAL PROFILE AND FVL+PT 3,03			RISK SCORE BASED ON THROMBO INCODE 7,52		
GENDER	H:1 M:2	2	GENDER	H:1 M:2	2
AGE	Years old	56	AGE	Years old	56
SMOKER		1	SMOKER		1
DIABETES		0	DIABETES		0
FAMILY HISTORY	NO: 0 YES: 1	1	FAMILY HISTORY	NO: 0 YES: 1	1
PREGNANCY		0	PREGNANCY		0
CONTRACEPTION		0	CONTRACEPTION		0
BMI	Kg/m ²	27	BMI	Kg/m ²	27
PROTHROMBIN (FII)	NO: 0 Hetero: 1	0	PROTHROMBIN (FII)		0
FV LEIDEN	Homo: 2	0	FV LEIDEN		0
			FV CAMBRIDGE		0
			FV HONG KONG		0
			FXII	NO: 0 Hetero: 1 Homo: 2	2
			FXIII		0
			SERPIN C1		0
			SERPIN A10		0
			A1 CARRIER		1

PATIENT CANDIDATE PROFILES




- Patients with a **family history** of VTE or hereditary thrombophilia.
- Patients with a **personal or family history** of VTE with **temporary conditions** that increase the risk of thrombosis.*
- Family members of an individual diagnosed with hereditary thrombophilia (**family study**).
- Patients being treated for venous thromboembolism to assess the **risk of re-thrombosis**.
- Patients with a VTE profile that suggests **hereditary thrombophilia**: VTE in patients under 45, recurring VTE, in unusual vascular areas, etc.
- **Women with issues of infertility** (repeated miscarriages/implantation failures to implant) with no identified cause or suspected thrombophilia.

* Such as: pregnancy, hormone-based contraception, hormone-replacement therapy, major surgery, paralysis or trips over 5 hours.


THROMBO INCODE REPORT

RISK SCORE AND MODULATION OF THE RISK

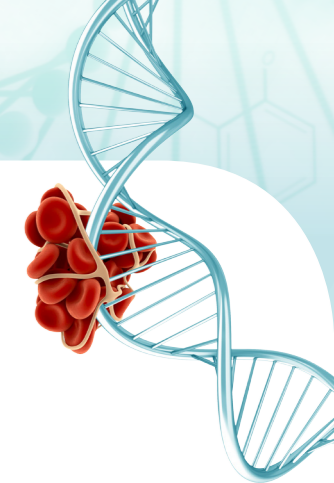

Assessment of risk for having venous thromboembolic events


According to the clinical and genetic information provided, your patient's current risk of having a thromboembolic event is 7.52.

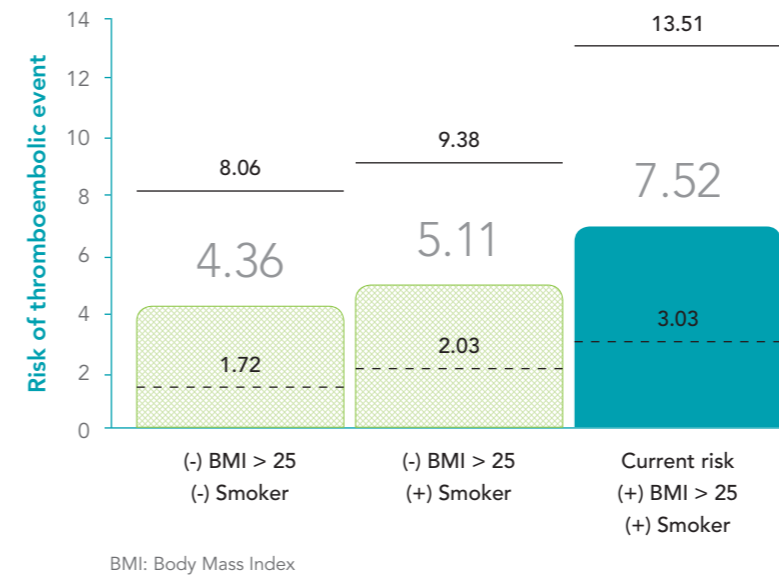


The patient's risk is compared to:

- The risk level of a similar patient without any additional genetic mutations (dotted line).
- The risk level of a patient with a similar clinical profile that is a heterozygote FVL carrier (solid line).

Evolution of patient's risk over time

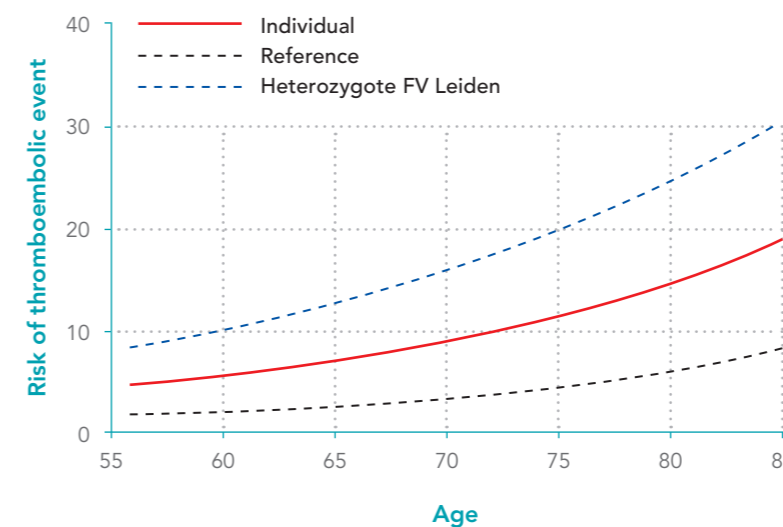

The patient's evolution of risk is compared to the same 2 theoretical profiles mentioned before.

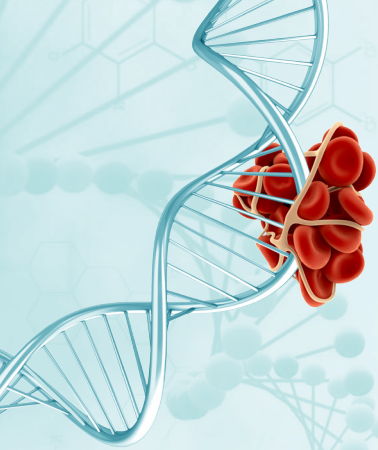

THROMBO inCode

DIAGNOSIS OF GENETIC THROMBOPHILIA PROFILE Greater sensitivity in detecting pro-thrombotic genetic mutations⁴

RISK SCORE FOR HAVING A VENOUS THROMBOEMBOLIC EVENT Overall risk assessment using both genetic and clinical data

PERSONALISED RECOMMENDATIONS REPORT Support in managing VTE risk of patients and their family members

COST/BENEFIT STUDY FOR VTE RISK ASSESSMENT Thrombo inCode is the predominant option (more effective and affordable than FVL and PT)²⁶

GEN inCode

Customer service
customerservice@genincode.com
(+34) 936 690 321
www.genincode.com

References: **1.** Bezemer ID et al. Predictive genetic variants for Venous Thrombosis: What's New?. *Semin Hematol* 44:85-92. **2.** Souto JC & Soria JM. Predicting individual risk of venous thrombosis. *Blood*. 2012 Jul 19; 120(3):500. **3.** Souto JC et al. Genetic susceptibility to thrombosis and its relationship to physiological risk factor: The GAIT study. *Am J Hum Genet*. 2000; 67:1452-1459. **4.** Macalalun P et al. Diagnosis and management of heritable thrombophilia. *BMJ*. 2014;4387:1-9. **5.** Soria JM et al. Multilocus Genetic Risk Scores for Venous Thromboembolism Risk. *J Am Heart Assoc*. 2014;3:e001060; originally published October 23, 2014. **6.** I.D. Bezemer et al. Predictive genetic variants. *Semin Hematol*. 44:85-92. **7.** Dölek B et al. Molecular analysis of factor V Leiden, factor V Hong Kong, factor II G20210A, methylenetetrahydrofolate reductase C677T, and A1298C mutations related to Turkish thrombosis patients. *Clin Appl Thromb Hemost*. 2007;13(4):435-8. **8.** Renner W et al. Prothrombin G20210A, factor V Leiden, and factor XIII Val34Leu: common mutations of blood coagulation factors and deep vein thrombosis in Austria. *Thromb Res*. 2000; 99(1):35-9. **9.** Young G et al. Influence of the factor II G20210A variant or the factor V G1691A mutation on symptomatic recurrent venous thromboembolism in children: an international multicenter cohort study. *J Thromb Haemost*. 2009 Jan;7(1):72-9. **10.** Bertina RM et al. Mutation in blood coagulation factor V associated with resistance to activated protein C. *Nature*. 1994; 369(6475):64-67. **11.** Williamson D et al. Factor V Cambridge: a new mutation (Arg306>Thr) associated with resistance to activated protein C. *Blood*. 1998; 91(4):1140-4. **12.** Chang WP et al. A novel mutation of Arg306 of factor V gene in Hong Kong Chinese. *Blood*. 1998; 91(4):1135-9. **13.** Emmerich J et al. Combined effect of factor V Leiden and prothrombin 20210A on the risk of venous thromboembolism—pooled analysis of 8 case-control studies including 2310 cases and 3204 controls. Study Group for Pooled-Analysis in Venous Thromboembolism. *Thromb Haemost*. 2001; 86(3):809-16. **14.** Poort SR et al. A common genetic variation in the 3'-untranslated region of the prothrombin gene is associated with elevated plasma prothrombin levels and an increase in venous thrombosis. *Blood* 1996;88:3698-3703 **15.** Souto JC et al. The prothrombin 20210A allele is the most prevalent genetic risk factor for venous thromboembolism in the Spanish population. *Thromb Haemost*. 1998; 80(3):366-9. **16.** Soria JM et al. A quantitative trait locus in the human factor XII gene influences both plasma factor XII levels and susceptibility to thrombotic disease. *Am J Hum Genet*. 2002; 70(3):567-574. **17.** Cochery-Nouvellon E et al. Homozygosity for the C46T polymorphism of the F12 gene is a risk factor for venous thrombosis during the first pregnancy. *J Thromb Haemost*. 2007; 5(4):700-7. **18.** Tirado I et al. Association after linkage analysis indicates that homozygosity for the 46C->T polymorphism in the F12 gene is a genetic risk factor for venous thrombosis. *Thromb Haemost*. 2004; 91:899-904. **19.** Morelli VM et al. ABO blood group genotypes and the risk of venous thrombosis: effect of factor V Leiden. *J Thromb Haemost*. 2005; 3(1):183-185. **20.** Ohira T et al. ABO blood group, other risk factors and incidence of venous thromboembolism: the Longitudinal Investigation of Thromboembolism Etiology (LITE). *J Thromb Haemost*. 2007; 5(7):1455-61. **21.** Tirado I et al. The ABO blood group genotype and factor VIII levels as independent risk factors for venous thromboembolism. *Thromb Haemost*. 2005; 93(3):468-74. **22.** Corral J et al. A nonsense polymorphism in the protein Z-dependent protease inhibitor increases the risk for venous thrombosis. *Blood*. 2006; 108(1):177-183. **33.** Tait RC et al. Prevalence of antithrombin deficiency in the healthy population. *Br J Haematol*. 1994 May;87(1):106-12. **23.** Corral J et al. Antithrombin Cambridge II (A384S): an underestimated genetic risk factor for venous thrombosis. *Blood*. 2007;109(10):4258-4263.401. **24.** Ariens RA et al. The factor XIII V34L polymorphism accelerates thrombin activation of factor XIII and affects cross-linked fibrin structure. *Blood*. 2000; 96(3):988-95. **25.** Catto AJ et al. Association of a common polymorphism in the factor XIII gene with venous thrombosis. *Blood*. 1999; 93(3):906-8. **26.** Rubio-Terrés C et al. Economic analysis of thrombo inCode, a clinical-genetic function for assessing the risk of venous thromboembolism. *Appl Health Econ Health Policy*. 2015 Apr;13(2):233-42

This material is only for use by healthcare professionals and is not meant for the general public. Distribution of this material in other countries shall be at the sole discretion and liability of official distributors of the products, who must be aware of and ensure compliance with current laws regarding distribution of materials in their country. GEN inCode shall not be held liable for any violation and/or failure to comply with local laws in other countries resulting from the actions of the distributors of the products.